



# ELASTICS

f-Ties™

Slingties™

Elastomerics

Intraoral Elastics



# Glide-Ties™

Glide-Ties™ reduce archwire binding, which may provide faster tooth movement and reduced treatment time. Elastomeric Glide-Ties™ are made with 'low friction technology' specifically targeted towards the reduction of frictional binding forces in orthodontic auxiliaries i.e., ligatures, separators, etc.

The super-tough FDA approved medical grade polyurethane is **latex-free**, odorless, tasteless, highly stain resistant and designed to deliver a uniform continuous force over a long period of time. They are injection molded in the conventional manner and then subjected to a hydro-thermal process that modifies the surface texture, creating a high degree of 'slipperiness'. This effect is permanent and significantly reduces treatment times.

ELASTICS

### Regular Glide-Ties™

24 Ties Per Stick  
1,008 Ties Per Pack  
24 color choices



### Mini Glide-Ties™

10 Ties Per Stick  
1,000 Ties Per Pack  
24 color choices



### Glide Chain Elastic

15 Foot Spool  
24 color choices



Continuous Short Long

### Regular Glide-Ties™ Kit

18 Popular Colors  
180 Sticks Per Kit  
4,320 Total Ties



Full Kit  
#300-200

\*Laminated Color Chart included with kit

Glide-Ties™ Dispenser  
#300-226



### Mini Glide-Ties™ Kit

24 Popular Colors  
480 Sticks Per Kit  
4,800 Total Ties



Full Kit  
#300-100+

\*Laminated Color Chart included with kit

Mini Glide-Ties™ Dispenser  
#300-126L



### Glide-Ties™ Combo Kit

12 Popular Colors  
120 Regular Sticks Per Kit  
12 Short Chains of matching color



Full Kit  
#300-300

\*Laminated Color Chart included with kit

Mini Glide-Ties™ Dispenser  
#700-062



Color	Regular 	Mini 	Chain  (continuous)	Chain  (short)	Chain  (long)
 Clear	#300-201	#300-101	#300-501	#300-301	#300-401
 Silver	#300-202	#300-102	#300-502	#300-302	#300-402
 White	#300-203	#300-103	#300-503	#300-303	#300-403
 Black	#300-204	#300-104	#300-504	#300-304	#300-404
 Red	#300-205	#300-105	#300-505	#300-305	#300-405
 Green	#300-206	#300-106	#300-506	#300-306	#300-406
 Royal Blue	#300-207	#300-107	#300-507	#300-307	#300-407
 Orange	#300-208	#300-108	#300-508	#300-308	#300-408
 Yellow	#300-209	#300-109	#300-509	#300-309	#300-409
 Purple	#300-210	#300-110	#300-510	#300-310	#300-410
 Periwinkle	#300-211	#300-111	#300-511	#300-311	#300-411
 Teal	#300-212	#300-112	#300-512	#300-312	#300-412
 Neon Green	#300-213	#300-113	#300-513	#300-313	#300-413
 Neon Pink	#300-214	#300-114	#300-514	#300-314	#300-414
 Aqua Blue	#300-215	#300-115	#300-515	#300-315	#300-415
 Light Pink	#300-216	#300-116	#300-516	#300-316	#300-416
 Light Green	#300-217	#300-117	#300-517	#300-317	#300-417
 Light Blue	#300-218	#300-118	#300-518	#300-318	#300-418
 Lilac	#300-219	#300-119	#300-519	#300-319	#300-419
 Midnight Blue	#300-220	#300-120	#300-520	#300-320	#300-420
 Turquoise	#300-221	#300-121	#300-521	#300-321	#300-421
 Burgundy	#300-222	#300-122	#300-522	#300-322	#300-422
 Dark Gray	#300-223	#300-123	#300-523	#300-323	#300-423
 Metallic Green	#300-224	#300-124	#300-524	#300-324	#300-424
 Pearl	#300-227	#300-127	#300-527	#300-327	#300-427
 Gold	#300-228	#300-128	#300-528	#300-328	#300-428
 Watermelon	#300-229	#300-129	#300-529	#300-329	#300-429
 Amethyst	#300-230	#300-130	#300-530	#300-330	#300-430
 Sparkle Blue	#300-231	#300-131	-	-	-
 Sparkle Green	#300-232	#300-132	-	-	-
 Sparkle Red	#300-233	#300-133	-	-	-
Assorted Colors	#300-225	#300-125	-	-	-

**Color Chart #000-005L**



**Chain Dispenser #700-065**



**Bulk Extraoral Elastics**

1/2" | 8 oz. | 1,000/pk | #350-608  
 1/2" | 16 oz. | 1,000/pk | #350-610



# Chroma Lastics

Orthodontic Elastics

ortho  
sense

## ADVANCED INTRA-ORAL ORTHODONTIC ELASTICS

Crafted in the USA, these elastics offer unparalleled quality and are available in a diverse range of options. Precision-engineered from premium medical-grade materials, ChromaLastics ensure reliable, uniform force for optimal orthodontic outcomes.

## ORTHODONTIC ELASTICS (LATEX + NON-LATEX)

All catalog numbers are given for latex. For non-latex, add "NL" to the catalog number when ordering

### PATIENT-SIZED QUANTITIES - 50 packets of 100 elastics

Made from the highest grade pure surgical latex tubing.

BAG COLOUR DIAMETER	PURPLE 1/8" (3.2mm)	GREEN 3/16" (4.7mm)	RED 1/4" (6.4mm)	BLUE 5/16" (7.9mm)	ORANGE 3/8" (9.5mm)	COLOUR
<b>2.5oz Light</b>	OS811-000	OS811-003	OS811-006	OS811-009	OS811-012	NATURAL
<b>4.5oz Medium</b>	OS811-001	OS811-004	OS811-007	OS811-010	OS811-013	NATURAL
<b>6.5oz Heavy</b>	OS811-002	OS811-005	OS811-008	OS811-011	OS811-014	NATURAL
<b>2.5oz Light</b>	OS811-175	OS811-176	OS811-177	OS811-178	OS811-179	NEON
<b>4.5oz Medium</b>	OS811-015	OS811-017	OS811-019	OS811-021	OS811-023	NEON
<b>6.5oz Heavy</b>	OS811-016	OS811-018	OS811-020	OS811-022	OS811-024	NEON

## BULK QUANTITIES - packets of 10,000

Standard elastics bagged bulk for special dispensing requirements.

DIAMETER	1/8" (3.2mm)	3/16" (4.7mm)	1/4" (6.4mm)	5/16" (7.9mm)	3/8" (9.5mm)	COLOUR
<b>2.5oz Light</b>	OS811-025	OS811-028	OS811-031	OS811-034	OS811-037	NATURAL
<b>4.5oz Medium</b>	OS811-026	OS811-029	OS811-032	OS811-035	OS811-038	NATURAL
<b>6.5oz Heavy</b>	OS811-027	OS811-030	OS811-033	OS811-036	OS811-039	NATURAL
<b>2.5oz Light</b>	OS811-170	OS811-171	OS811-172	OS811-173	OS811-174	NEON
<b>4.5oz Medium</b>	OS811-040	OS811-042	OS811-044	OS811-046	OS811-048	NEON
<b>6.5oz Heavy</b>	OS811-041	OS811-043	OS811-045	OS811-047	OS811-049	NEON

## BULK QUANTITIES - packets of 100,000

Standard elastics bagged bulk for special dispensing requirements.


DIAMETER	1/8" (3.2mm)	3/16" (4.7mm)	1/4" (6.4mm)	5/16" (7.9mm)	3/8" (9.5mm)	COLOUR
<b>2.5oz Light</b>	OS811-200	OS811-203	OS811-206	OS811-209	OS811-212	NATURAL
<b>4.5oz Medium</b>	OS811-201	OS811-204	OS811-207	OS811-210	OS811-213	NATURAL
<b>6.5oz Heavy</b>	OS811-202	OS811-205	OS811-208	OS811-211	OS811-214	NATURAL
<b>2.5oz Light</b>	OS811-230	OS811-231	OS811-232	OS811-233	OS811-234	NEON
<b>4.5oz Medium</b>	OS811-220	OS811-222	OS811-224	OS811-226	OS811-228	NEON
<b>6.5oz Heavy</b>	OS811-221	OS811-223	OS811-225	OS811-227	OS811-229	NEON

## EXTRA-ORAL ELASTICS

DIAMETER	1/4"	5/16"	3/8"	1/2"	COLOUR
<b>8oz</b>	OS811-200	OS811-203	OS811-206	OS811-209	NATURAL
<b>16oz</b>	OS811-201	OS811-204	OS811-207	OS811-210	NATURAL

## Slingties™

Ortho Arch Slingties™ have many innovative features. They are designed for moving the teeth to the wire instead of moving the wire to the teeth, which negates the reciprocal forces on the adjacent teeth. They are best used in cases that require bringing down high cuspids, moving teeth into space, closing diastemas, midline corrections and anchoring teeth together. The 6mm and 9mm Slingties may be stretched up to three times their size with 100% return for up to eight weeks. The translucent color resist staining with aesthetic or stainless steel brackets.

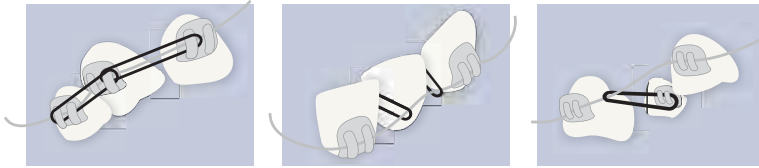

Package contains:  
 (104) 6mm Slingties™  
 (104) 9mm Slingties™  
**#300-000**



## Super Chain

This is a proprietary specialty polyurethane material that has superior memory and very high elastic recovery as compared to conventional polyurethane materials (CPM) used widely in orthodontics.



**Silver**  
 15 ft. Spool  
 Short #300-001  
 Long #300-002  
 Cont. #300-003

**Clear**  
 15 ft. Spool  
 Short #300-004  
 Long #300-005  
 Cont. #300-006



## Elastomeric Separators

These separators are radiopaque and made from elastomeric material. They are designed to separate teeth more rapidly due to it's faster recovery and are easily placed with the Deluxe Separating Pliers.

1,000 pk  
**#325-302**

Deluxe Separating Pliers  
**#800-010**




## Safety Knob Separators

This uniquely designed separator has two raised knobs to prevent it from drifting sub-gingival. It is easily stretched and placed between two opposing teeth, providing great bounce back properties to create the desired space for your band.

6 Stars (480 Modules Light Blue Thin)  
**#325-300**


6 Stars (480 Modules Purple)  
**#325-301**

Deluxe Separating Pliers  
**#800-010**

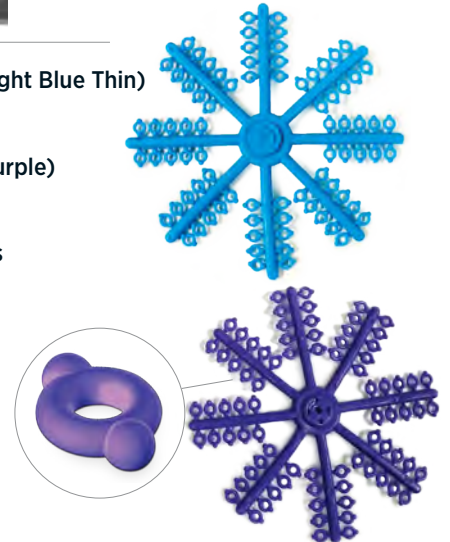

## Elastomeric Tissue Guard

This strong elastic tubing will fit over archwires or small springs to prevent irritation and accommodates wire up to .025 diameter

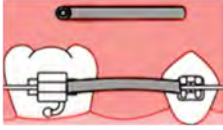

Gray		Clear	
10 ft. Spool		10 ft. Spool	
.018	#340-309	.018	#340-304
.027	#340-310	.027	#340-305
.045	#340-311	.045	#340-306



## Elastomeric Thread & Tubing


This thread holds knots securely. Solid tubing offers greater force.

Tube - Gray		Tube - Clear	
25 ft. Spool		25 ft. Spool	
.025	#340-210	.025	#340-205
.030	#340-230	.030	#340-225
Thread - Gray		Thread - Clear	
25 ft. Spool		25 ft. Spool	
.025	#340-110	.025	#340-105
.030	#340-130	.030	#340-125



## Rotation Wedges

These rotation wedges are designed to lock in to place adjacent to the tie wing of bracket, providing a fulcrum between the wire and bracket for rotation.

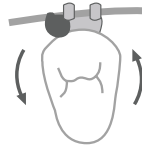

Gray	100/pk	#270-022
Clear	100/pk	#270-021



## Elastic Attacher

This handy tool is an aid to assist patients in the placement and removal of elastics.



Assorted Colors	100/pk	#800-026
-----------------	--------	----------



## Mathieu Needle Holder w/ Safety Latch - Ultra Fine Tip

Best Mathieu for Ligatures



#800-121

Ultra-Fine Tip



Safety Latch